


MicroPeak

The precision Peak Flow
Meter with an immaculate
pedigree

Micro
Medical



With EN13826
European Scale

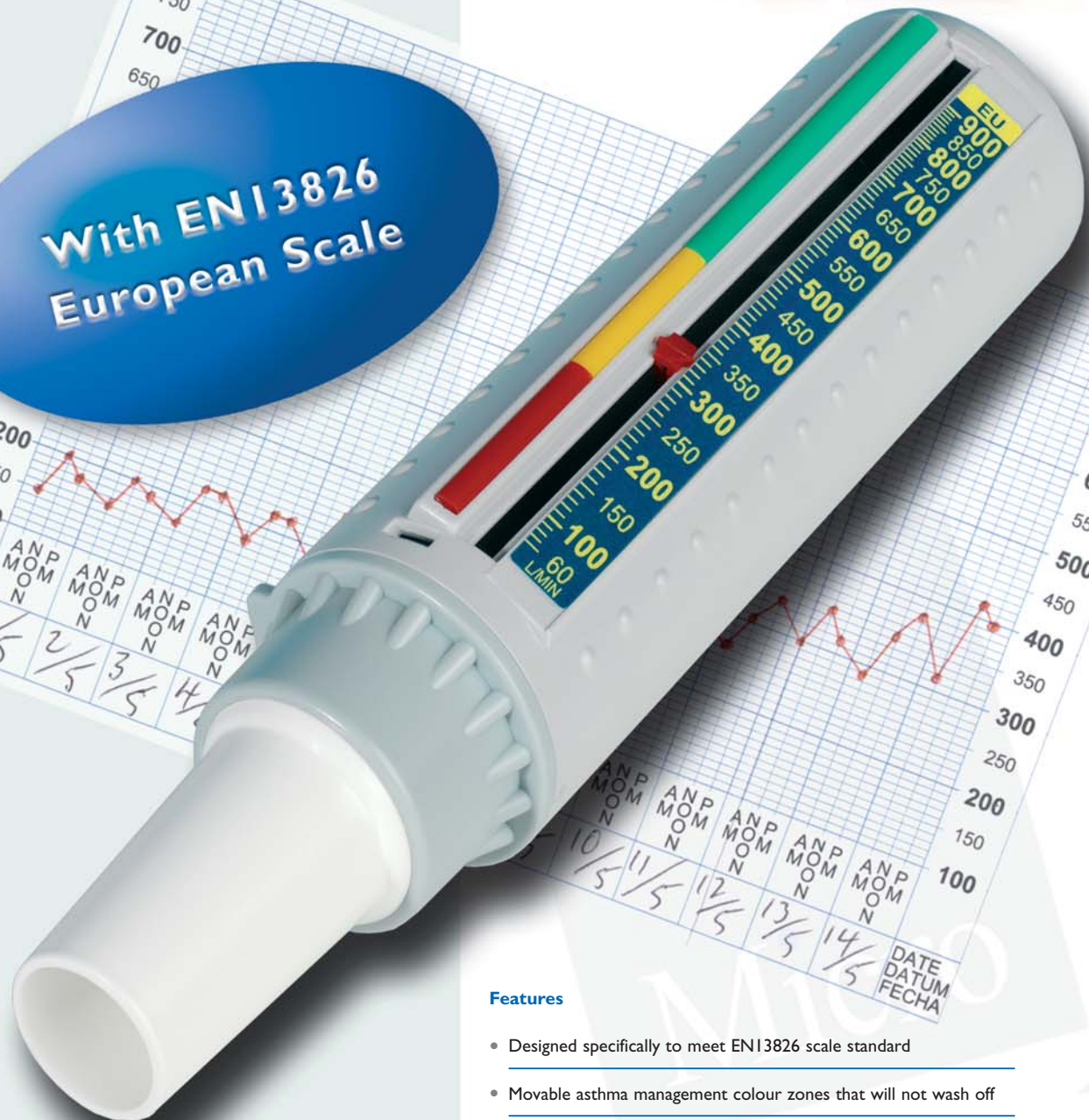

Features

- Designed specifically to meet EN13826 scale standard
- Movable asthma management colour zones that will not wash off
- Large, high visibility scale
- In-built internal check valve, avoiding cross infection concerns
- Meets and exceeds all international performance standards
- Low cost
- Customisation possible for promotional purposes
- Available with either European or ATS scale

25 years experience of design, testing and improving the world's leading Peak Flow Meters (PFM) has led to the new generation MicroPeak.


Designed to meet the new European (EN13826) standards in performance and quality including a high visibility scale and in-built

Colour-Zone asthma management system that won't wash off (see illustration below).

Unlike other PFM's the advanced design of MicroPeak means dismantling is avoided, thus protecting the integrity of the internal mechanism, also an internal check valve avoids cross infection concerns.

MicroPeak comes complete with instructions and a 4 week recording chart, and can be offered fully customised for promotional purposes too.

When accuracy is important in clinical asthma management, then compromise is not an option, choose MicroPeak the PFM with an immaculate pedigree.

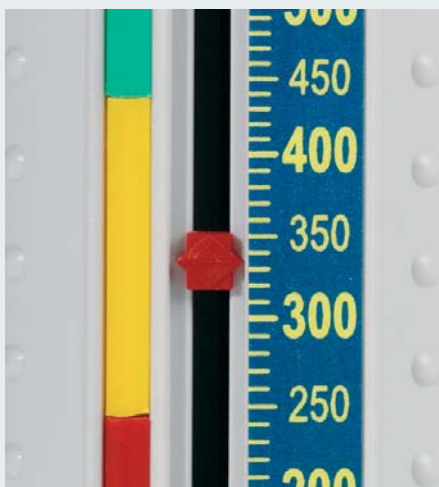

(Colour-Zone asthma management system)


MicroPeak Cat. No. MPE8200EU is available on UK drug tariff

MicroPeak EU scale ▶

MicroPeak ATS scale ▼

Technical Specification

- Range: 60-900 l/min (EN13826, (Cat. No. MPE8200EU))
60-900 l/min (ATS, (Cat. No. MPE7200))
- Scale increment 10 l/min
- Accuracy \pm 5% or 10 l/min
- Repeatability \pm 2%
- Calibration 100% (individual unit)
- Frequency response 100Hz @ 350 l/min
- Dynamic state critically damped @ 350 l/min
- Standard compliance ERS statement 1997 (Eur Respir J 1997:9 Suppl.24)
AS/NZS 4237:1994 NAEP Guidelines 1991
- Internal check value expiration only

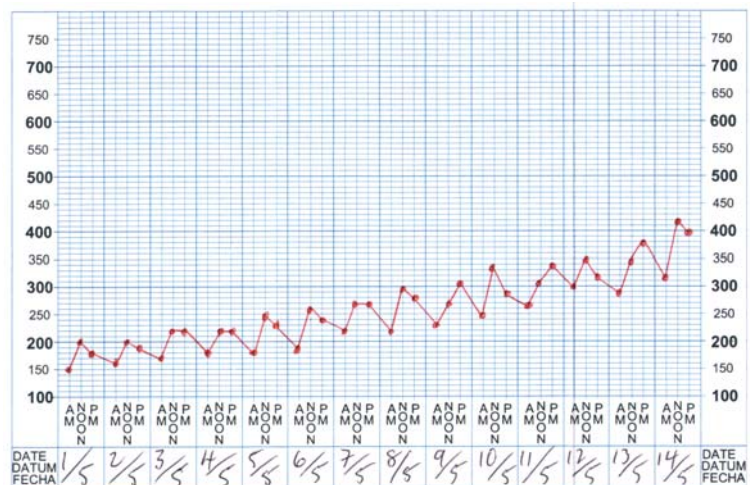

Peak flow record chart showing recovery from acute asthma with pronounced morning dips.

Micro Medical Limited

PO Box 6, Rochester,
Kent, ME1 2AZ, UK

Telephone 01634 893500

Fax 01634 893600

International +44 1634 893500

Email sales@micromedical.co.uk

www.micromedical.co.uk



These Peak Flow Meters are part of an extensive range of respiratory measurement devices manufactured by Micro Medical Ltd and are offered as Catalogue No MPE8200EU MicroPeak European scale or Catalogue No MPE7200 MicroPeak ATS (American Thoracic Society) scale.

Micro Medical Ltd pursues a policy of continuing improvement in design, production and performance of its products, the right is therefore reserved to vary details at any time without notice.