



Pulse Oximetry

The world's #1 oximeter

9550 Onyx II



All-in-one
fingertip
pulse oximeter

lightweight

reliable

Vital Signs Within Reach.™



Product Highlights

- ✓ Simple - Automatically turns on for fast, easy spot-checks
- ✓ Superior Accuracy - Trusted by clinicians worldwide
- ✓ Efficient - New turn off mechanism allows for 2500 spot-checks or 21 hours of continuous use
- ✓ Cost Effective - No additional sensors needed
- ✓ Enhanced Design - Enclosed spring; smaller; wider viewing angle; easy battery access

Nonin, the world leader in fingertip oximetry, has designed years of clinical experience into the Onyx II. The electronics and sensor are incorporated into the world's smallest fingertip oximeter, providing a cost effective solution for spot-checks and short term monitoring. Onyx II – not just a new look. Superior performance and unparalleled reliability put vital signs within reach anywhere, anytime.

✓ Where People Make the Difference



Vital Signs Within Reach

The Onyx II is ideal for use in any situation where a fast accurate reading of blood oxygen saturation and pulse rate is needed. Never search for a pulse oximeter or sensor again.



Proven Performance/ Simple Operation




✓ **Small, Lightweight** - Carry the Onyx II anywhere – on a lanyard, in a carrying case or in your pocket. It weighs only 2 oz.

✓ **Large LED display** - The widest viewing angle available for optimum visibility of patient readings, day or night.



✓ **Battery Door** - Enhanced design for convenient battery access.

✓ **Pulse Quality Indicator** - Brilliant, easy-to-see LED provides a quick and accurate assessment of patient readings.

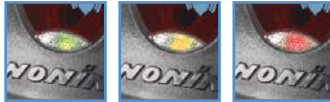

✓ **Energy Efficient** - A new, turn-off mechanism allows two AAA batteries to provide approximately 2500 spot checks or 21 hours of continuous operation.



✓ **Pediatric to Adult Finger Sizes** - Accommodates a wide range of finger thicknesses, 0.3-1.0 inch (8 mm to 26 mm).

✓ **Durable** - So rugged and reliable, it can be used in the most demanding environments.

Everything You Need



Standard Package Items



Optional Carrying Cases

Please contact your authorized distributor for more information.

Specifications

Oxygen Saturation Range (%SpO ₂)	0% to 100%
Pulse Rate Range	18 to 321 pulses per minute
Displays	Two 3-digit LEDs, 7 segment LED (light emitting diodes)
Numeric Displays	Tricolor LED segments
Pulse Strength	
Measurement Wavelengths	
Red	660 nanometers
Infrared	910 nanometers
Accuracy	
Blood Oxygen Saturation (%SpO ₂ ±1SD) ^a	70 - 100% ± 2 digits
Pulse Rate	20 - 250 bpm ±3%
Temperature	
Operating	+32° to +122°F (+0° to +40°C)
Storage/Transportation	-22° to +122°F (-30° to +50°C)
Humidity	
Operating	10% to 90% noncondensing
Storage/Transportation	10% to 95% noncondensing
Altitude	
Operating Altitude	Up to 40,000 feet
Hyperbaric Pressure	Up to 4 atmospheres
Power Requirements	Two 1.5 volt AAA-size batteries
Battery Life	
Operating	Approximately 21 hours of continuous operation ^b
Storage	4 years, minimum ^c
Dimensions	1.31 x 2.22 x 1.27 in. (3.32 x 5.63 x 3.22 cm)
Weight	1.89 oz (53.5 g) with batteries
Warranty	2 years

Classifications per IEC 60601-1

Degree of Protection	Type BF-Applied Part
----------------------	----------------------

^aS.D. (Standard Deviation) is a statistical measure; up to 32% of the readings may fall outside these limits.

^bUsing two AAA-size alkaline batteries; based on approximately 2500 spot checks calculated at 30 seconds per spot check.

^cUsing two AAA-size alkaline batteries.

Quality systems registered to ISO 9001:1994/EN 46001:1996, ISO 13485:1996, and CE marking per Annex II, Clause 3 of EC Directive No. 93/42/EEC concerning medical devices.



Specifications are subject to change without notice.

Your Partner of Choice

With over 18 years of experience and dedication to the design and support of noninvasive monitoring devices, Nonin has helped many medical professionals meet clinical and economic objectives. Trusted by clinicians worldwide, countless tests are performed every day on thousands of Nonin pulse oximeters in more than 125 countries. Nonin's dedication to technological leadership, precision manufacturing and uncompromised customer support ensures quality products and service you expect.



www.nonin.com

Nonin Medical, Inc.
13700 1st Avenue North
Plymouth, MN 55441-5443, USA
Phone: +1.763.553.9968 Toll Free: 1.800.356.8874
Fax: +1.763.553.7807
info@nonin.com